


# CONFIGURING A PHYSIOTEL DIGITAL IMPLANT FOR BPR ANALYSIS IN PONEMAH V5.20

The following steps explain how to enable a soft channel for a PhysioTel™ Digital implant in order for it to be used with Ponemah's BPR analysis module to determine respiration from blood pressure.

1. With the CLC connected to the system, launch **Ponemah v5.20**.
2. Select **Hardware | Hardware | Edit Configuration...**

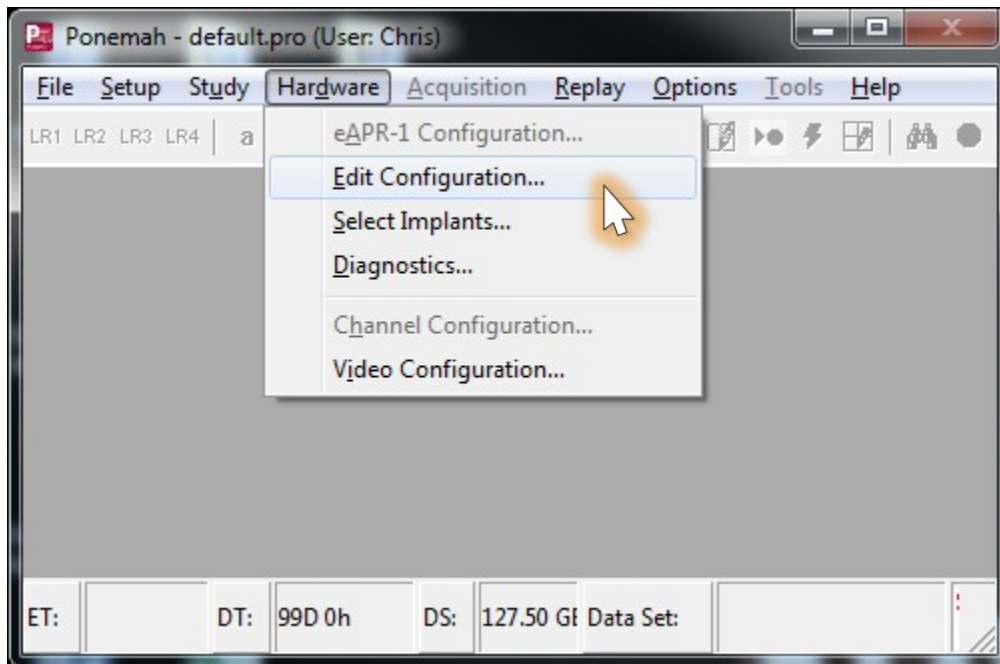

3. Configure the PhysioTel Digital implants using the *Configure Sources* wizard.
4. From the list view on the left, expand the implant of interest to display its channels.
5. Select **Signal Strength** from the list of channels.  
*Note:* Signal Strength will be configured for use with BPR and will no longer be able to be used to collect the Signal Strength waveform.
6. Click the checkbox associated with **Enabled**.
7. Change the *Description* field to **RespPress**.

**Configuration**

Ready Configuration

- DSICLC-671063
  - 673565
    - Pressure 1
    - Pressure 2
    - Biopotential
    - Temperature
    - Accelerometer X
    - Accelerometer Y
    - Accelerometer Z
    - Battery Life Remaining
    - Activity
    - RespPress**
    - Battery Days Remaining
  - 660732
  - 673568
  - 669304

**Soft Channel Details**

Description: RespPress

Implant Name: 673565

Channel Source: SignalStrength

Enabled: ☒

8. **Repeat** this process for each implant that will use the BPR analysis module.
9. Click **Save & Exit** when done.
10. Configure the protocol as normal.

Channel Input Setup

| Input                       | Analysis | Label    | Units   |
|-----------------------------|----------|----------|---------|
| 1 - (673565 - Pressure 1)   | LVP      | LVP      | mmHg    |
| 2 - (673565 - Pressure 2)   | BP       | Pressure | mmHg    |
| 3 - (673565 - Biopotential) | ECG      | ECG      | mV      |
| 4 - (673565 - Temperature)  | TEMP     | Temp     | Celsius |
| 5 - (673565 - Activity)     | ACT      | Activity | Counts  |
| 6 - (673565 - RespPress)    | BPR      | RespPres | mmHg    |
| 7 - (660732 - Pressure 1)   | BPR      | VP       | mmHg    |
| 8 - (660732 - Pressure 2)   | CBF      | ressure  | mmHg    |
| 9 - (660732 - Biopotential) | CVOL     | CG       | mV      |
| 10 - (660732 - Temperature) | CYS      | emp      | Celsius |
| 11 - (660732 - Activity)    | DIS      | activity | Counts  |
|                             | ECG      |          |         |
|                             | ECGRO    |          |         |
|                             | EMG      |          |         |

An extra channel should appear in Ponemah's **Channel Input Setup** for each implant. Configure this channel to use the **BPR** analysis module, and configure the *Labels* and *Units* accordingly.

**Note:** If using **Auto Configure Protocol** function to automatically define protocol settings, Ponemah will setup the Respiration channel as **Signal Strength** using the **RAW** analysis module. Simply change the analysis module to **BPR**, then adjust the *Labels* and *Units* accordingly.