


CONFIGURING A PHYSIOTEL DIGITAL IMPLANT FOR BPR ANALYSIS IN PONEMAH V5.20

The following steps explain how to enable a soft channel for a PhysioTel™ Digital implant in order for it to be used with Ponemah's BPR analysis module to determine respiration from blood pressure.

1. With the CLC connected to the system, launch **Ponemah v5.20**.
2. Select **Hardware | Hardware | Edit Configuration...**

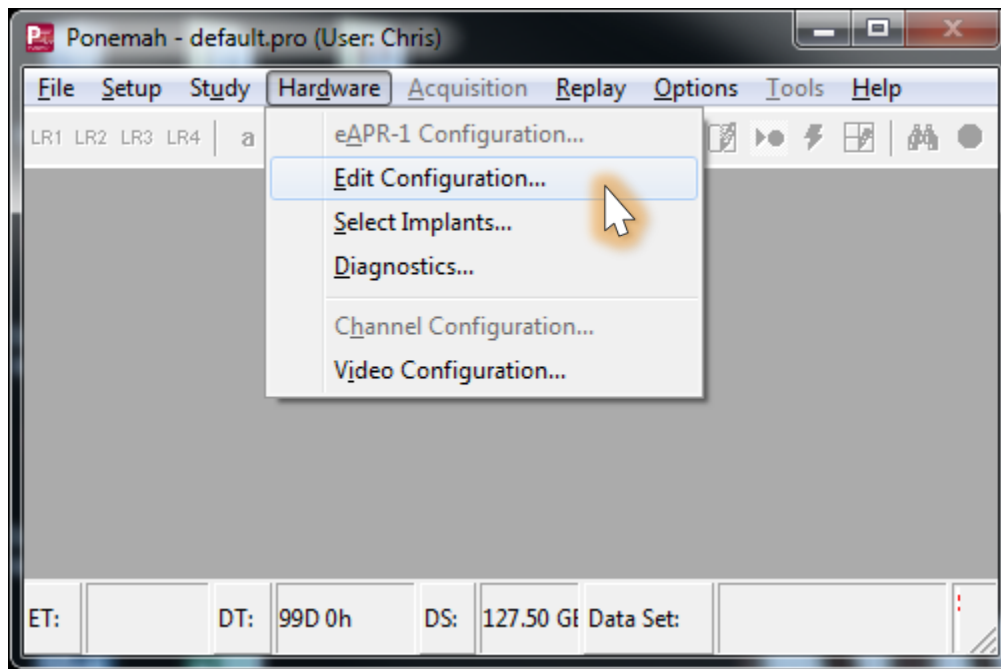

3. Configure the PhysioTel Digital implants using the *Configure Sources* wizard.
4. From the list view on the left, expand the implant of interest to display its channels.
5. Select **Signal Strength** from the list of channels.

Note: Signal Strength will be configured for use with BPR and will no longer be able to be used to collect the Signal Strength waveform.

6. Click the checkbox associated with **Enabled**.
7. Change the *Description* field to **RespPress**.

Configuration	Soft Channel Details
Ready Configuration ● DSICLC-671063 ▲ ● 673565 Pressure 1 Pressure 2 Biopotential Temperature Accelerometer X Accelerometer Y Accelerometer Z Battery Life Remaining Activity RespPress Battery Days Remaining ▶ ● 660732 ▶ ● 673568 ▶ ● 669304	Description: RespPress Implant Name: 673565 Channel Source: SignalStrength Enabled: ☒

8. **Repeat** this process for each implant that will use the BPR analysis module.
9. Click **Save & Exit** when done.

Configure the protocol as normal.

Channel Input Setup				
	Input	Analysis	Label	Units
up	1 - (673565 - Pressure 1)	LVP	LVP	mmHg
	2 - (673565 - Pressure 2)	BP	Pressure	mmHg
up	3 - (673565 - Biopotential)	ECG	ECG	mV
	4 - (673565 - Temperature)	TEMP	Temp	Celsius
Protocol Head	5 - (673565 - Activity)	ACT	Activity	Counts
Setup	6 - (673565 - RespPress)	BPR	RespPres	mmHg
Setup	7 - (660732 - Pressure 1)	BPR	VP	mmHg
	8 - (660732 - Pressure 2)	CBF	ressure	mmHg
	9 - (660732 - Biopotential)	CVOL		
	10 - (660732 - Temperature)	CYS	CG	mV
rt	11 - (660732 - Activity)	DIS	emp	Celsius
on		ECG	ctivity	Counts
		ECGRO		
		EMG		

An extra channel should appear in Ponemah's **Channel Input Setup** for each implant. Configure this channel to use the **BPR** analysis module, and configure the *Labels* and *Units* accordingly.

Note: If using **Auto Configure Protocol** function to automatically define protocol settings, Ponemah will setup the Respiration channel as **Signal Strength** using the **RAW** analysis module. Simply change the analysis module to **BPR**, then adjust the *Labels* and *Units* accordingly.

BLACK-BINOPTIC electronic tap

Mains supply, H. 170mm

Spout with smooth interior

Reference: **388035**

Specification description

Deck-mounted electronic basin tap:

Independent IP65 electronic control unit.

230/12V mains supply with transformer.

Flow rate pre-set at 3 lpm at 3 bar, can be adjusted from 1.4 to 6 lpm.

Scale-resistant flow straightener.

Adjustable duty flush (pre-set flush ~60 seconds every 24 hours after last use).

Active, infrared presence detection, sensor at the end of the spout optimises detection.

Chrome-plated brass body.

Matte black chrome finish.

PEX flexible with filter and solenoid valve M3/8".

Fixing reinforced by 2 stainless steel rods.

Anti-blocking security.

Spout with smooth interior and low volume of water (reduces bacterial development).

Drop height: 170mm for semi-recessed washbasins.

Suitable for people with reduced mobility.

10. 30-year warranty.